

PRODUCT NO.: 7080



Dispenser, MD Pen Airless 6.5mL, Capacity 0.10mL/click

Description: The new portfolio of MD Dispensers was engineered for accuracy with a professional structure, yet designed to convey a cosmetic appeal. The MD Dispenser will maintain the integrity of a formulation, allowing users to dose preparations easily and accurately. The first of its kind on the market, this dispenser features a UV resistant view window for added convenience.

Features:

- Airless dispensing
- Slim, sleek pen-type look
- Silicone tip for smooth application
- Metered dispensed volume
- Excellent evacuation
- Minimal priming required
- UV-resistant
- Protection from oxidation
- Hygienic use, no bacterial contamination
- Latex-, BPA- and PVC-free

Technical Specifications:

Capacity (mL)	6.5
Dispensed Volume (mL)	0.10

USA | Canada | Australia | International
medisca.com


Please note that the information contained herein is for informational purposes only. The information is based on recognized compendia, manufacturers' instructions and/or other authoritative sources and texts. MEDISCA takes no responsibility for the validity or accuracy of this information. Customer MUST refer to USP, FDA, Health Canada, TGA and/or their State Board to ensure they comply with the specific regulations and ensure that appropriate procedures are in place.

Precision (% RSD)*	< 5.0
Material of Construction	Piston: Polyethylene Inner Bottle: Polypropylene Outer Bottle: Polypropylene Cap: Polypropylene Nozzle: Polypropylene Tip: Silicone
Dimensions, L x Ø (cm / in)	13.9 x 2.0 / 5.5 x 0.8
Color	White with view window
Filling Method	3 piece assembly: 1. Top-Filled Pen 2. Snap-Cap Spring Pump 3. Outer Casing
UV-Resistance	Bottle: Yes View Window: Yes

* Value established using Medisca cream bases.

Filling Instructions:

MD Pen

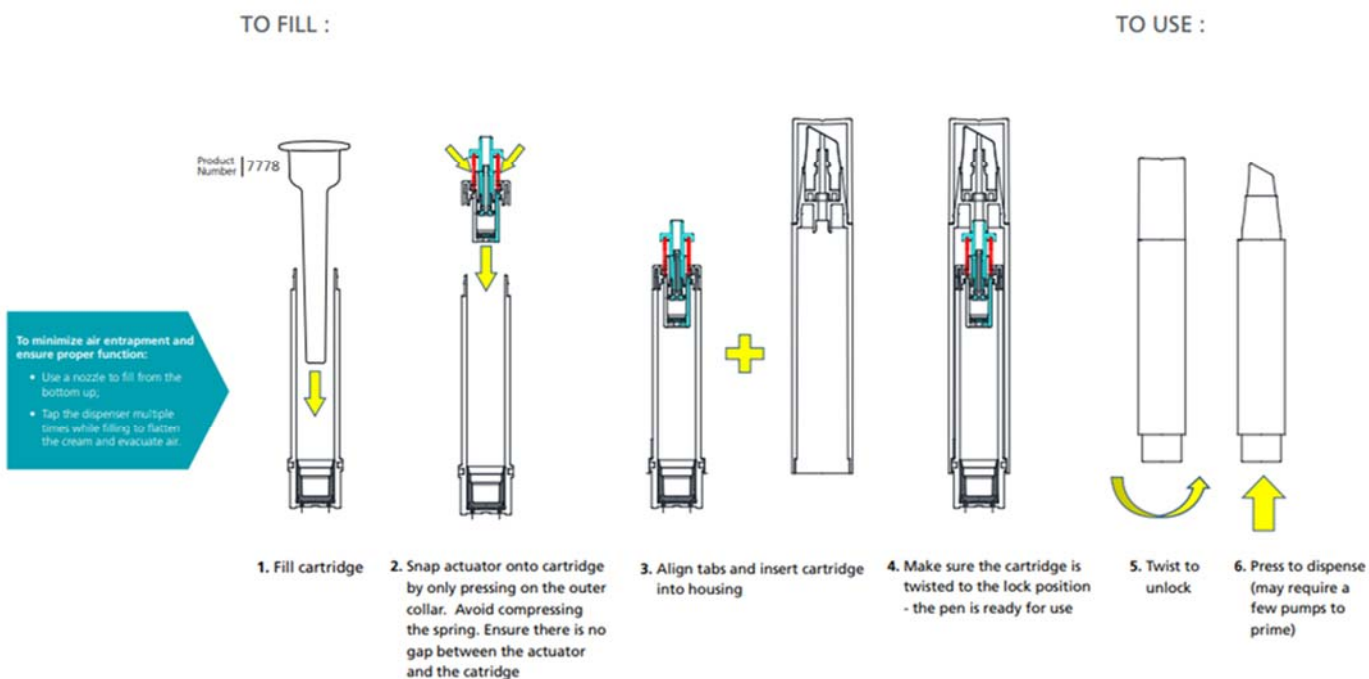

USA | CANADA | AUSTRALIA PACIFIC | EUROPE
WWW.MEDISCA.COM

Please note that the information contained herein is for informational purposes only. The information is based on recognized compendia, manufacturers' instructions and/or other authoritative sources and texts. MEDISCA takes no responsibility for the validity or accuracy of this information. Customer MUST refer to USP, FDA, Health Canada, TGA and/or their State Board to ensure they comply with the specific regulations and ensure that appropriate procedures are in place.