








PRODUCT NO.: 5622

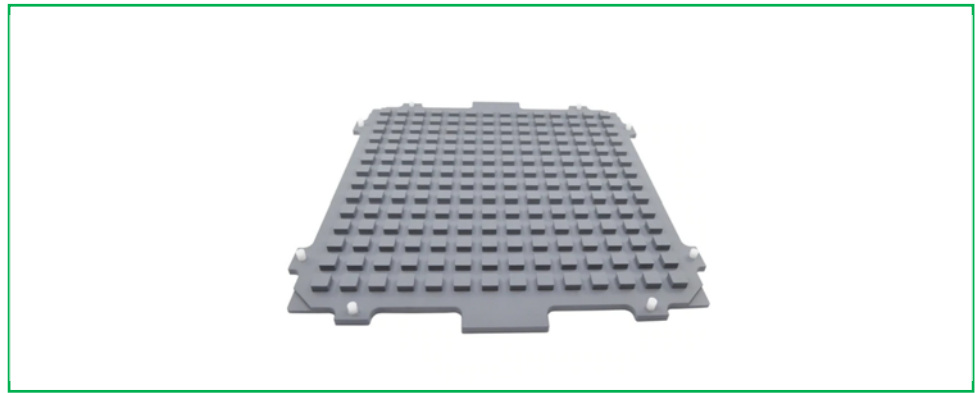

Ejektor Plate, Medi X-TABB

Description: The Ejektor Plate is specifically designed for easy and efficient dose removal from the Cavity Plate in Medi X-TABB System. Simply place Ejektor Plate on a flat solid surface with bosses facing up, line up the cavity plate between the guides and press down on the cavity plate to eject the doses. It is constructed of durable and lightweight aluminum that is coated with the latest US FDA and EU compliant non-stick technology that provides an easy-to-use and easy-to-clean surface.



FEATURES:

-  Made of durable and lightweight aluminum coated with US FDA and EU compliant non-stick technology
-  Easy-to-use and easy-to-clean surface
-  Hard coating resists scratches and reasonable wear when used in accordance with instruction provided
-  Teflon guides for easy and consistent placement

TECHNICAL SPECIFICATIONS:

Dimensions (inches) (Approx.)	13 x 13
Material	Hard coated Aluminum
Total Number of bosses	196

Additional Information

-  Used with Cavity Plates (5623 and 5628) and Quarter Score Plate (5627)
-  Sold separately. Not included with Medi X-TABB Mold (RDT System) (Code: 5650)

USA | CANADA | AUSTRALIA PACIFIC | EUROPE
WWW.MEDISCA.COM

USA | CANADA | AUSTRALIA PACIFIC | EUROPE
WWW.MEDISCA.COM

Please note that the information contained herein is for informational purposes only. The information is based on recognized compendia, manufacturers' instructions and/or other authoritative sources and texts. MEDISCA takes no responsibility for the validity or accuracy of this information. Customer MUST refer to USP, FDA, Health Canada, TGA and/or their State Board to ensure they comply with the specific regulations and ensure that appropriate procedures are in place.